

Adjustable OA Defiance®

Custom Osteoarthritis Knee Brace

The adjustable OA Defiance® is designed for the active OA patient that does not want to compromise their high level of activity. This lightweight custom brace provides durable support for moderate to severe levels of osteoarthritis. The bilateral hinge bars offer support for ligament instabilities and work in conjunction with a patented telescoping condyle pad which puts the load control in the hands of the patient.



A1




Through an exclusive arrangement with dj Orthopedics, VQ OrthoCareSM offers the **VQ BY DJ** brand of orthopedic products. These products add the value of VQ OrthoCare's exclusive sales, service and patient care network to the DonJoySM line of quality products.



Telescoping Pad



Key

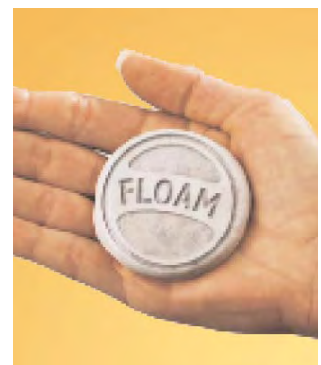

Floam™

Patient controlled telescoping condyle pad allows patient to fine-tune load according to activity level

Adjustable OA Defiance®

Custom Osteoarthritis Knee Brace

Product Profile

Features and Benefits

The Adjustable OA Defiance® brace is designed to stabilize the joint and alter the knee's mechanics with its three-point loading system. A force is applied across the center of the knee joint and a reaction force responds at the cuffs in the opposite direction, redistributing the weight away from the painful compartment.

- Base load correction of 1-5 degrees built into the frame
- Patient controlled telescoping condyle pad allows patient to fine-tune load according to activity level
- Lightweight, low-profile, four-walled carbon composite frame molded to accommodate unique varus/valgus angulations
- Floam™ condyle pad with 100% memory foam for increased patient comfort and compliance
- Patented Four-Points-of-Leverage™ system provides excellent ligament stability
- Custom-made utilizing the DonJoy CCMi measurement system or cast mold
- Suede frame and strap pads standard

OPTIONS

- Available in standard or short lengths
- ACL, PCL or CL versions
- 40 custom color or graphic designs available



Ordering Information

Each Adjustable OA Defiance® is custom made using measurements taken with the CCMi Mark II or casted in weight bearing

Size	Thigh	Calf
Minimum	12 1/2" (31.8 cm)	10 1/4" (26 cm)
Maximum	38 1/8" (96.5 cm)	26 1/4" (66.4 cm)
Accessories		
3-1107-0-00000	Floam™ condyle kit	Universal
3-0177-0-00000	Chamois liner	Custom
3-0757-X-06000	Chamois strap pads	S, M, L

Indications

Mild to severe unicompartmental osteoarthritis

Warning

Knee Bracing

If you experience any pain, swelling, sensation changes or any unusual reactions while using this product, contact Patient Care at VQ OrthoCare, 1-800-452-7993 or consult a physician.

Neoprene supports should not be worn by individuals with known susceptibility to dermatitis or allergies to neoprene byproducts.

CAUTION

Federal law restricts this device to sale by or on the order of a practitioner licensed by the law of the State in which he/she practices to use or order the use of the device.

Notice

While every effort has been made in state-of-the-art techniques to obtain the maximum compatibility of function, strength, durability and comfort, there is no guarantee that injury will be prevented through the use of this product.



18011 Mitchell South, Irvine, CA 92614
 Call 800.266.6969 Fax 800.821.8012
 www.vqorthocare.com
 customerservice@vqorthocare.com

Floam is a registered trademark of Mattel, Inc.
 Four-Points-of-Leverage is a trademark and
 Defiance and DonJoy are registered trademarks of dj Orthopedics, LLC
 VQ OrthoCare is a service mark of VisionQuest Industries, Inc. © 2005