

# Saunders Cervical HomeTrac<sup>®</sup>



## Cervical Traction

Our sleek new design features self-adjusting neck wedges that eliminate the need for awkward adjustment knobs. The new wedges provide closer, more secure contact for a perfect fit every time. Lie down, fasten the head strap and let the HomeTrac's comfortable wedges self-adjust by rotating on their base. Never has home cervical traction been simpler.

- Ensures optimum angle for traction treatment
- Provides therapeutic force up to 50 pounds
- Replicates clinical traction
- Self-adjusting neck wedges
- Allows traction in supine position
- Ensures continuity between home and clinical treatments

# Saunders Cervical HomeTrac®

## Cervical Traction

### Product Profile

#### Features and Benefits


- Directs traction forces toward the occiput (back of the head) thereby preventing compression of the TMJ through the chin
- Applies up to 50 pounds of traction
- Ensures safety - user has total control of force at all times
- Enables easy home traction without assistance
- Provides a comfortable and cost-effective option to continuing clinical traction treatments
- Includes instructional video, user's guide and deluxe carrying case
- Patented pump design incorporates an easy-to-use "click-in-place" mechanism for pumping, sustaining and releasing the traction



- Gauge reads in pounds/kg of traction to provide feedback to the user and to promote consistency in treatment



- Self-adjusting neck wedges rotate through a 20° arc and are removable for easy cleaning or replacement

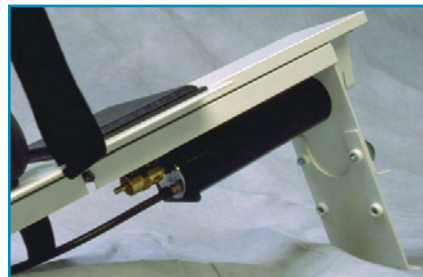

- Adjustable extension foot allows for additional traction angles of 20° and 25°



Before



During

- The cervical spine before traction and during traction applied by the Saunders Cervical HomeTrac®

#### Indications

- Pain in back or neck
- Herniated or protruding discs
- Degenerative disc disease
- Posterior facet syndrome
- Sciatica
- Joint hypomobility
- Sub acute joint inflammation
- Nerve root impingement
- Paraspinal muscle spasm

#### Contraindications

- Acute or traumatic injury
- Spinal instability
- Fractures
- Extruded disc fragmentation
- Severe rheumatoid arthritis
- Spinal cord compression
- Infections or inflammatory disease
- Serious pathology, disease or malignancy
- Spreading and/or aggravation of patient's symptoms

#### Code

#### Description

CTRAC DLX	Saunders Cervical HomeTrac
100625-R2	Replacement neck wedges

#### Precautions

It is recommended that the force be increased gradually over several sessions. Slight, post-treatment muscle soreness in the neck is common; however excessive soreness or an increase in peripheral symptoms may result if the force is increased too quickly.

#### CAUTION

Federal law restricts this device to sale by or on the order of a practitioner licensed by the law of the State in which he/she practices to use or order the use of the device.



18011 Mitchell South, Irvine, CA 92614  
800.266.6969 Fax 800.821.8012  
www.vqorthocare.com  
customerservice@vqorthocare.com

Made in the U.S. Patent Pending.  
HomeTrac is a registered trademark of The Saunders Group Inc.  
VQ OrthoCare is a service mark of VisionQuest Industries, Inc. © 2005